



MACH7 Integration Manager Quick Start Guide

Jose Sanchez - 2024-12-19 - Comments (0) - MACH 7

Integration Manager Quick Start Guide

You may receive a prompt that there is a new integration manager version available for download when launching images in DX Studio via the integration manager where already installed on a local client. You will follow same steps as listed in below snippet for downloading and updating the integration manager.

*Steps to be followed a MACH7 VNA upgrade.

*Please clear your local Edge browser cookies cache before accessing the Mach7 apps.

Integration Manager Quick Start Guide

The Mach7 Integration Manager is a downloadable application which enables integration between the Universal Worklist and the eUnity viewer. Settings such as preferred browser and saving layouts can be configured in the application.


Workstation Requirements

- All workstations running Integration Manager will require .NET version 4.5 or higher be installed.

Installation

To download and install the M7T Integration Manager:

1. From within Diagnostic Studio, select **Download M7T Integration Manager** from the Tools dropdown menu on the Main Worklist screen.

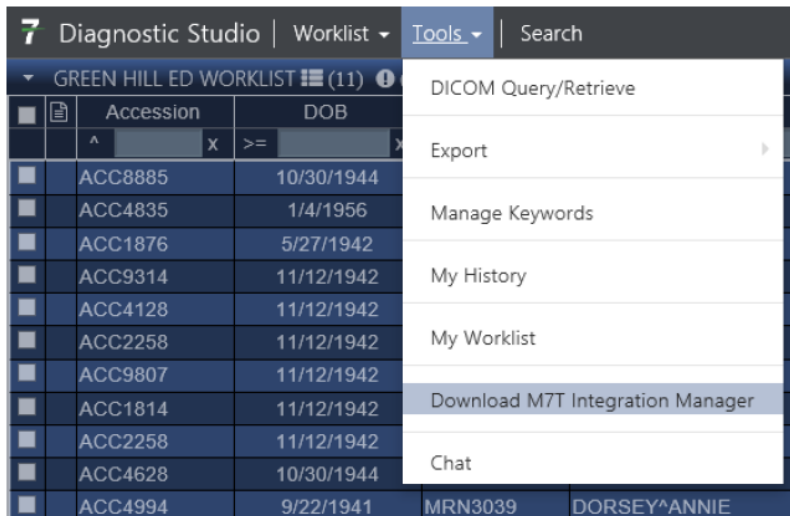

2. Choose whether to run or save and launch the M7T Integration Manager installer when prompted.
3. Run the Mach7 Integration Manager Setup will launch and continue through the installation wizard to install the Mach7 Integration manager. When prompted, enter the IP or Fully Qualified Domain Name (FQDN) of the server hosting the Mach7 Enterprise Imaging Platform.

Tags
Integration Manager
MACH7

Attachments

- [MACH7 - Integration Manager Quick Start Guide.docx \(178.58 KB\)](#)
- [Integration Manager Quick Start Guide.pdf \(190.38 KB\)](#)